

APOGEE SPECTRORADIOMETERS | PS-100, PS-200, and PS-300

Absolute spectral measurement across a wide wavelength range

Product Specifications

Features

Three Wavelength Options 350 to 1000 nm, 300 to 830 nm, or 300 to 1000 nm.

Complete Package

Includes spectroradiometer, two meter fiber-optic cable, cosine-corrected detector, USB cable, USB drive with required drivers and software (compatible with all 32-bit and 64-bit Windows operating systems), and shoulder bag (functions as a carrying case and field measurement pack). A reflectance probe and reflectance standard are available as accessories.

Portable Field Measurements


Features a small and lightweight design with a rugged housing and no moving parts. Spectroradiometer is powered through the USB port on a computer allowing mobile field measurements.

Typical Applications


Measurement of spectral output of radiation sources for plant or human lighting, reflectance and transmittance measurements of natural and synthetic surfaces and materials (often plant leaves and canopies), and absorptance measurements of chemical samples.

	PS-100	PS-200	PS-300	
Irradiance Calibration Range	350 to 1000 nm	300 to 850 nm	300 to 1000 nm	
Wavelength Sensitivity	350 to 1150 nm	190 to 850 nm	220 to 1100 nm	
Wavelength Resolution	1 nm	0.85 nm	1.5 nm	
Detector Type	CCD, 2048 pixel			
Grating Type	Holographic & Ruled, 600 g/nm	Holographic and aberration-corrected, 590 g/nm		
Digitizer	16-bit			
Signal to Noise Ratio	1000:1			
Stray Light	0.1 % at 435 nm, 0.5 % at 600 nm	0.02 % at 435 nm, 0.2 % at 200 nm	0.02 % at 435 nm, 0.2 % at 220 nm	
Measurement Repeatability	Less than 1 %			
Irradiance Calibration Uncertainty	± 10 %			
Detector Integration (Exposure) Range	1 ms to 65 s			
Directional (Cosine) Response	± 5 % at 80° zenith angle			
Software	Windows compatible, included			
Computer Interface	USB 2.0			
Power Requirement	100 mA at 5 V DC, supplied via USB code			
Operating Temperature	0 to 60 C			
Optical Cable	2 m armored fiber-optic			
Base Unit Size	25 mm x 75 mm x 125 mm	I 60 mm v 100 mm v 150 mm		
Mass	500 g	900 g		
Warranty	1 year against defects in materials and workmanship			



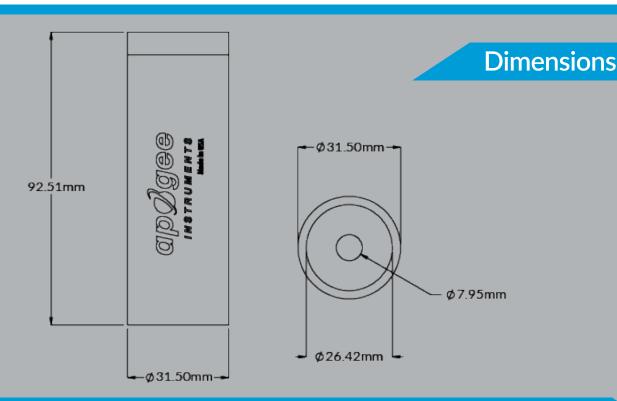

Screenshot of Measurement Made in SpectraWiz Software

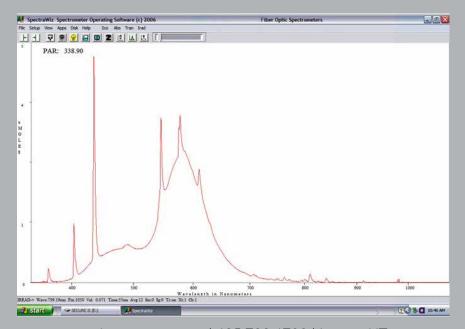

The detector collects light for the spectroradiometer via fiber-optic cable. The spectroradiometer then transmits its measurements to a PC where custom software (included) displays the spectrum of the light source. Live data can be saved as a data file suitable for work in other programs, such as Excel or SigmaPlot.